


TANITA®

Monitoring Your Health

Dual Frequency Body Composition Monitor with integrated Bluetooth technology

innerScan®
DUAL

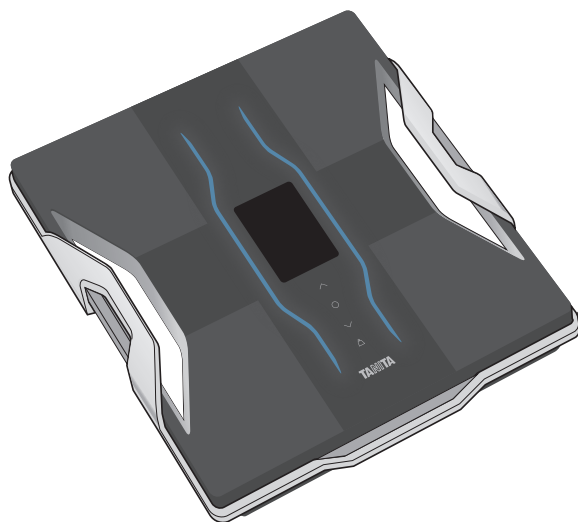

Model: RD-901



Instruction Manual

Please read this Instruction Manual carefully and keep for future reference.



Table of Contents

Introduction 3

Features and Functions 3

For Your Safety 4

Getting Accurate Readings 5

Preparations before Use 6

Connecting to an iPhone 7

Setting (Resetting) Date and Time 8

Setting and Storing Personal Data in the Memory 9

Taking Measurements 11

Measurement Results 13

Specifying the Personal Data Number for Measurement 14

Programming the Guest Mode 15

Measuring Weight Only 15

Body Composition Guide 16

Troubleshooting 22

Specifications 23

Introduction

Thank you for purchasing this Tanita Body Composition Monitor.

This monitor is one of a wide range of home health care products produced by Tanita.

This Instruction Manual will guide you through the setup procedures and outline the key features of

this monitor. Additional information on healthy living can be found on our website www.tanita.asia

Tanita products incorporate the latest clinical research and technological innovations.

All data is monitored by the Tanita Medical Advisory Board to ensure accuracy.



Safety Precautions

Do not use the Body Composition Monitor if you have an electronic medical implant such as a pacemaker, as it passes a low-level electrical signal through the body which may interfere with the operation of the implant.

Pregnant women should only use the weight function. Other functions are not intended for use when pregnant.

Important Notes for Users

This Body Composition Monitor is intended for adults and children (age 5-17) who have inactive to moderately active lifestyles, and adults with athletic lifestyles.

Thanks to major investment in the latest BIA Technology and sports science research,

TANITA has upgraded the Athlete Mode function to make measurements more accurate and suitable for a wider range of users.

Who should use Athlete Mode?

Adults aged 18 years and over who either:

- Train or exercise for 12 hours or more a week and have been doing so for at least six months.
- Are body builders.
- Are professional athletes who want to monitor their progress at home.
- Have a lifetime history of fitness and who used to do more than 12 hours a week but do less now.

Recorded data may be lost if the unit is used incorrectly or is exposed to electrical power surges.



Tanita cannot accept any responsibility for losses incurred due to the loss of recorded data.

Tanita cannot accept any responsibility for damages or losses incurred through the use of this product or any third-party claims.


Note: Body fat percentage estimates will vary with the amount of water in the body, and can be affected by dehydration or over hydration due to such factors as alcohol consumption, menstruation, illness, intense exercise, etc.

Features and Functions


Front

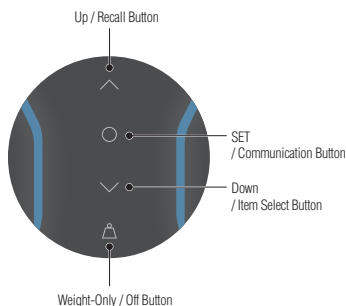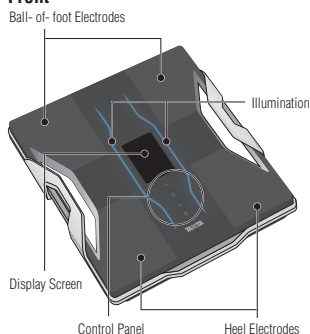

Back

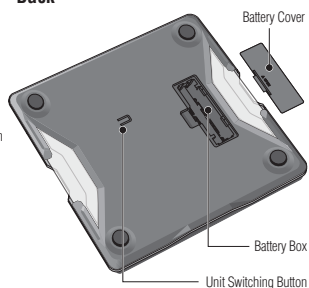


Quick Start Guide

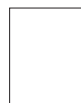

Accessories

AA-Size Batteries (4)

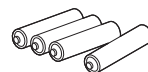

For Your Safety

This section explains precautionary measures to be taken to avoid injury to the users of this device and others, and to prevent damage to property. Please familiarise yourself with this information to ensure safe operation of this equipment.

	Warning	Failure to follow instructions highlighted with this mark could result in death or severe injury.
	Caution	Failure to follow instructions highlighted with this mark could result in injury or damage to property.
		This mark indicates actions that are prohibited.
		This mark indicates instructions that must always be followed.

Warning

 Prohibited	This monitor must not be used by people with pacemakers or other medical implants. This monitor passes a weak electrical current through the body which could interfere with and cause the malfunction of electrical medical implants, resulting in serious harm to the user.
	Interpretation of measurements (e.g. evaluation of measurements and formulation of exercise programmes based on the measurements) must be performed by a professional. Implementing weight loss measures and exercise programmes based on self-analysis could be harmful to your health. Always follow the advice of a qualified professional.
	Do not step on the edge of the platform. Do not jump onto the monitor. Doing so may cause you to fall or slip, resulting in serious injury.
	Do not place this monitor on a slippery surface such as a wet floor. Doing so may cause you to fall or slip, resulting in serious injury. It may also cause water to get inside the monitor resulting in damage to internal components.
	Do not drop any objects onto the monitor. Avoid excessive impact or vibration to the monitor. Doing so may damage the glass, resulting in injury.

Caution

 Prohibited	Do not stand on the display screen or control panel.
	Do not insert your fingers into gaps or holes. Doing so may result in injury.
	Do not use rechargeable batteries. Do not use old batteries together with new batteries, or a mix of different types of batteries at the same time. Doing so may cause battery fluid to leak or the batteries to become excessively hot and rupture, resulting in damage to the monitor or injury.
	Never allow an infant or small child to carry the monitor. The child may drop the monitor, resulting in injury.
	This monitor is intended for home use only. This monitor is not intended for professional use including hospitals, medical or fitness facilities. It is not equipped with the components required for heavy usage. Using the monitor in a professional location will invalidate the warranty.

Place the monitor on a hard, flat surface where there is minimal vibration to ensure safe and accurate measurement.

Avoid excessive impact or vibration to the monitor.

Doing so may damage the monitor or cause it to malfunction, or may cause setting to be erased.

Do not use or store this monitor in the following locations.

- Locations that are exposed to direct sunlight, near to heaters, or exposed directly to air from air conditioning units, etc.
- Locations where there are severe temperature changes or a large amount of vibration.
- Locations with high levels of humidity or moisture.

Doing so may damage the monitor or cause it to malfunction.

Do not step on the monitor when your body is wet (immediately after a shower or bath, etc.)

Be sure to thoroughly dry your body and feet before using the monitor. Failure to do so may cause you to slip or fall, resulting in injury. Also, the monitor cannot take accurate measurements if your body or feet are wet.

Keep away from water.

Do not use boiling water, benzene, thinner, etc. to clean the monitor.

If the monitor becomes dirty, soak a soft cloth in water or neutral household detergent, wring it out thoroughly and wipe the monitor clean, then use a dry cloth to wipe the monitor dry.

Do not attempt to disassemble the monitor. Do not modify this equipment in any way.

Do not wash the monitor in water.

Do not submerge this monitor in water.

Doing so may damage the monitor or cause it to malfunction.

Be sure to insert the batteries with the polarities (+, -) in the correct position.

If the batteries are incorrectly inserted, battery fluid may leak and damage the floor. (If you do not intend to use this monitor for a long period of time, it is advisable to remove the batteries before storage.)

Do not use the monitor near other products such as microwave ovens that emit electromagnetic waves.


Doing so may cause the monitor to malfunction, preventing accurate measurements.

Another person should assist persons with disabilities who may not be able to take a measurement alone.



Always...

Getting Accurate Readings

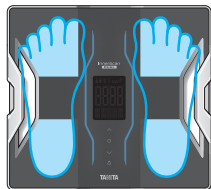

To ensure accuracy, readings should be taken without clothing and under consistent conditions of hydration. Always remove your socks or stockings, and be sure the soles of your feet are clean before stepping on the measuring platform. Be sure that the arch of each foot is aligned with the center of the measuring platform. Don't worry if your feet appear too large for the unit, accurate readings can still be obtained if your toes overhang the platform. It is best to take readings at the same time of day. Try to wait about two hours after getting up, eating, or strenuous exercise before taking measurements. While readings taken under other conditions may not have the same absolute values, they are accurate for determining the percentage change as long as the readings are taken in a consistent manner. To monitor your progress, compare your weight and body fat percentage readings taken under the same conditions over a period of time.


Note: An accurate reading will not be possible if the soles of your feet are not clean, if your knees are bent or if you are sitting down.

Preparations Before Use

Inserting the Batteries

Note:


- Be sure that the batteries are inserted in the correct orientation. If the batteries are inserted incorrectly, fluid may leak from the batteries and damage the floor.
- If you do not intend to use this unit for a long period of time, remove the batteries before storage.
- The included batteries from factory may have decreased energy levels over time.

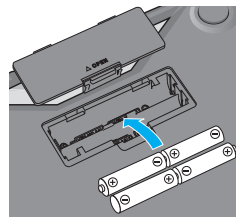

Set your location code

Different locations across the world have different gravitational pulls. These differences affect the accuracy of the weight measurement. You can obtain the most accurate weight readings by selecting the correct gravity setting on the scale according to your geographical location.

1. When the power is off, press and hold the Δ button.
2. Refer to the map enclosed, and press the $\wedge \vee$ buttons to select your location code.
3. Press the $\bigcirc$ button to enter and save the setting.

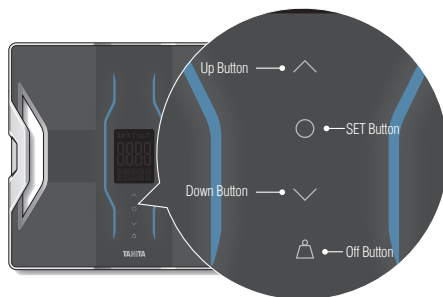

Changing the Weighing Unit

Note: Only for use with models that have the unit switching button.

The weighing unit can be changed using the switch on the back of the monitor.

Press the Δ button when the display is blank. When "0.00kg" (or "0.0lb", "0.0stlb") appears on the display, press the Unit Switching button underside of the monitor to switch the weighing unit.


Note: If weighing unit is set to pounds or stone-pounds, the height programming mode is automatically set to feet and inches. Similarly, if kilograms is selected, height is automatically set to centimeters.

Connecting to an iPhone

1. Turn the Bluetooth setting on, and use the URL to download the Health Planet app.
<http://www.tanita.asia/innerscandual/rd901.htm>

Note:

- See the Tanita website for details on supported models.
- If you do not wish to use the dedicated app, configure the settings manually on the body composition monitor. (P.8)
- You need to register your iPhone again if you replace your iPhone.

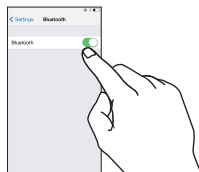


2. Launch the app. You need to register as a member to use the Health Planet app. After registering, login with your registered details.



3. Select "Input Data", then select "Input from BCM".

Note:

- When taking measurements while holding your iPhone, enter the weight that you want to subtract into the "Smartphone's weight" section in the dedicated app before taking measurements.
- The weight that you entered to be subtracted from your weight is displayed.
- The "Smartphone's weight" can be entered up to -2.00kg (-4.4lb, -4lb).
- The maximum weight (weighing capacity) is reduced by the amount entered for the "Smartphone weight".

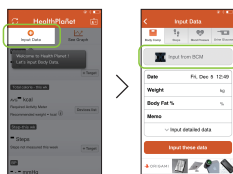

4. Check the Body Composition Monitor is off, and tap the "Add Device" button.

Note:


- The approximate communication distance is 5m in an unobstructed location.




5. Follow the instructions on the app, and press and hold the  (communication) button on the Body Composition Monitor for 3 seconds.

Note:

- When registering your body composition, place the Body Composition Monitor on a hard, flat floor and step on it with bare feet.

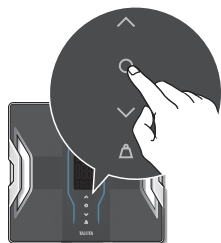

6. Follow the instructions on the app and measure your body composition. When measurement is complete, the connection to the iPhone is disconnected.

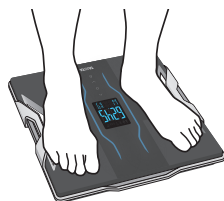


The screen displays and instructions may differ depending on the app specifications.

Setting (Resetting) Date and Time (Not required if connection to the iPhone is complete)

1. Set Date

While the power is off, press the $\bigcirc$ button to start setting.
Press the $\wedge \vee$ buttons to set year and then press the $\bigcirc$ button.
Set the month and date in the same way.

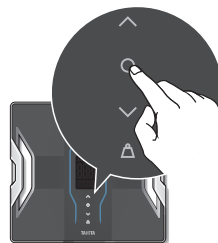

Note:

- When the power is switched on for the first time, the date and time setting screen can also be displayed by pressing the $\wedge \vee$ buttons.
- When you have replaced the batteries, the date and time setting screen can also be displayed by pressing the $\wedge \vee$ buttons.



2. Set Time

Press the $\wedge \vee$ buttons to set hour and then press the $\bigcirc$ button.
Set the minutes in the same way.
The set date and time are displayed in order, and then the power automatically turns off.



Note:

- If the Body Composition Monitor is not operated for 60 seconds or more, the power will turn off automatically. Settings will not be saved in this case.
- If you make a mistake or wish to turn the monitor off before you complete the settings, press the $\triangle$ button to turn the monitor off. Settings will not be saved in this case.
- To change the date or time settings, press and hold the $\vee$ button for 3 seconds when the power is off, and then follow the procedure above to change the date or time settings.

Setting and Storing Personal Data in the Memory (Not required if connection to the iPhone is complete)

The Body Composition Monitor can only be operated if data has been programmed into one of the personal data memories. Only weight measurements can be taken if you use the Body Composition Monitor without registering personal data.

1. Select a Personal Data Number

Press the $\bigcirc$ button to turn on the monitor. Press the $\wedge \vee$ buttons to select a personal data number (1, 2, 3 or 4).

Press the $\bigcirc$ button to confirm the personal data number.

Note:

- The monitor turns off automatically if it is not operated for 60 seconds.
- If you make a mistake or want to turn the monitor off before you have finished programming it, press the $\triangle$ button to force quit. Partially registered data is saved on the monitor.
- The backlight color changes according to the personal data number.

1:Blue

2:Pink

3:White

4:Light blue



2. Set Your Birthday

Set the date of your birthday (Year / Month / Day) using the $\wedge \vee$ buttons.

Press the $\bigcirc$ button to enter each item.

Note:

Year settings can be configured from 1900 onwards.



3. Select Female or Male

Use the $\wedge \vee$ buttons to scroll through Male [♂], female [♀], male athlete [♂🏃] and female athlete [♀🏃] settings, then press the $\bigcirc$ button.

Who should use Athlete Mode?

Adults aged 18 years and over who either:

- Train or exercise for 12 hours or more a week and have been doing so for at least six months.
- Are body builders.
- Are professional athletes who want to monitor their progress at home.
- Have a lifetime history of fitness and who used to do more than 12 hours a week but do less now.



4. Set Your Height

The default height is 170cm (5' 7.0") (height range is 90cm - 220cm or 3'0.0" - 7'3.0").

Use the $\wedge \vee$ buttons to set your height and then press the $\bigcirc$ button.

The display shows "0.00kg" (or 0.0lb, 0.0stlb) and the monitor is ready for measurement.



Setting and Storing Personal Data in the Memory

5. Measuring Your Body Composition

Check "0.00kg" (0.0lb, 0.0stlb) is shown on the display, then step on the monitor.
Registered personal data and the measurement results are displayed in order, then the power switches off automatically.

Note:

If you step onto the monitor before 0.00kg is displayed, "Err,-OL" is displayed and the power turns off. Set the personal information again.

Also, if you do not step onto the monitor within 60 seconds of "0.00kg"(0.0lb, 0.0stlb) appearing, the power turns off automatically. In this case, you need to set and store personal data again.

Note:

The monitor needs to be connected to the iPhone in order to use the dedicated app.

Connect to the iPhone, and register the Body Composition Monitor as shown in the instructions on P.7.



If you change the birthday and/or gender, "CLr" is displayed on the monitor.



Select "YES" using the $\wedge \vee$ buttons and then press the $\bigcirc$ button, to delete the existing personal data enabling you to register new data.



If you select "no" using the $\wedge \vee$ buttons and then press $\bigcirc$ button, no new data is saved (existing data is left unchanged) and power turns off.



Taking Measurements

1. When the monitor power is off, select "Input Data" in the app and then follow the instructions in the app.
When "0.00kg" (0.0lb, 0.0stlb) is displayed, step onto the monitor.



2. Measurement data is automatically sent to the iPhone.
If there is unsent data stored on the monitor, that data is also sent to the iPhone.


Note:

- When taking measurements while holding your iPhone, enter the weight that you want to subtract into the "Smartphone's weight" section in the dedicated app before taking measurements.
- The weight that you entered to be subtracted from your weight is displayed.
- The "Smartphone's weight" can be entered up to -2.00kg (-4.4lb, -4lb).
- The maximum weight (weighing capacity) is reduced by the amount entered for the "Smartphone weight".



If you do not have your iPhone when taking measurements

Simply step onto the monitor. The auto-recognition feature will automatically recall your personal data. Your monitor will store up to 10 sets of measurements.
Saved data is sent to the iPhone together with new measurements the next time the iPhone is connected to the monitor when taking measurements.

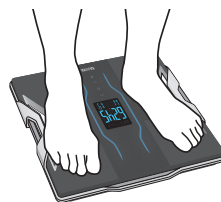

Taking Measurements

You can start taking measurements when you have finished entering your personal data.

1. Step onto the Monitor

Step onto the monitor when the power is off. The monitor starts measuring automatically.



2. View Your Results

After taking measurements, your results are automatically displayed in the order listed on the next page.

After the last result data is displayed, the weight is displayed again and then the monitor automatically turns off.

Note:

- Do not step off until “0000” disappears.
- The backlight color changes according to the recognized personal data number. (Blue/Pink/White/Light blue)
- When the monitor is connected to the iPhone and measurement data is saved in the monitor,  lights in the personal data display and the measurement result display screen.
- If there are more than 6 sets of data stored, the number of saved sets of data is displayed after the personal data display lights.
- Up to 10 sets of measurement data can be stored in the monitor.
- If you have set a nickname in the dedicated app, the nickname is displayed instead of the personal data number when personal data is displayed.

If the wrong personal number is displayed

Press the  button to select your personal number while readings or personal data is being displayed.

Then press the  button to display your own personal data and measurement results.

Situations where mis-recognition may occur

- When measuring people who have similar weights and body shapes
- When measurements have not been taken for a few days and there is a reasonable change in body weight
- When measuring directly after having the monitor stood upright, carrying the monitor or storing it with items placed on top

If the monitor continues to mis-recognize, specify your personal data number manually before measurement. (P.14)

Measurement Results

Weight



Body Mass Index



Body Fat %



Muscle Mass



Physique Rating



Bone Mass



Visceral Fat



Basal Metabolic Rate



Metabolic Age



Total Body Water %



An assessment of your measurement results are indicated by the color of the backlight.

Assessment Color(Applicable age)	●	●	●
BMI (used when Weight/BMI are displayed) (18-99)	Obese	Over / Under	Healthy
Body fat percentage (5-99)	Obese	Over / Under	Healthy(-) / Healthy(+)
Visceral fat level (18-99)	High	Slightly High	Average
Muscle mass (18-99)	Under	-	Average / High
Basal Metabolic Rate (used when BMR/Metabolic Age are displayed) (18-99)	Under	-	Average / High
Bone Mass (18-99)	Under	Slightly Under	Average / High

Note:

- For children (age 5-17), the monitor only displays weight, BMI, body fat percentage, muscle mass and bone mass.
- For children (age 4 or less) and adults (age over 100), the monitor only displays the weight.
- "Err" is displayed if the body fat percentage exceeds the measurable range.
- You can use the ∇ button to scroll through the results.
- Press and hold the ∇ button for fast scrolling.
- After reading the results that you want to see, release the button to see the rest of the results.
- The monitor automatically turns off after 40 seconds.
- Press the $\triangle$ button to cancel measurement or turn the power off.
- Do not move the monitor for approximately 3 seconds after taking measurements and stepping off.
To prevent accidents such as falling off the monitor, always step off the monitor before operating the buttons.

Specifying the Personal Data Number for Measurement

Follow these instructions if the monitor continues to mis-recognize your number.

1. When the power is off, press the $\wedge$ button to start measuring.
Select your personal number using the $\wedge/\vee$ buttons and then press the $\bigcirc$ button.



2. Step onto the monitor when "0.00kg" (0.0lb, 0.0stlb) appears in the display.
The monitor turns off automatically after displaying the measurement results.



Note:

- When the monitor is connected to the iPhone and measurement data is saved in the monitor,  lights in the personal data display screen and the measurement result display screen.
- If there are more than 6 sets of data stored, the number of saved sets of data is displayed after the personal data display lights.
- Up to 10 sets of measurement data can be stored in the monitor.
- If you step onto the monitor before 0.00kg (0.0lb, 0.0stlb) is displayed, "Err,-OL" is displayed and the power turns off during measurement.

Programming the Guest Mode

Guest mode allows you to program the monitor for one-time use without using a personal data number.


Press the $\wedge$ button to select "GUEST", then press the $\bigcirc$ button.

Next, follow the steps below:

Set the age, gender (male/ female/ male athlete/ female athlete) and height.

Step onto the monitor after "0.00kg" (0.0lb, 0.0stlb) is displayed.

The readings are displayed for 40 seconds.

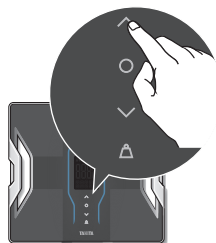

Note:

The age range is from 5 to 99.

Who should use Athlete Mode?

Adults aged 18 years and over who either:


- Train or exercise for 12 hours or more a week and have been doing so for at least six months.
- Are body builders.
- Are professional athletes who want to monitor their progress at home.
- Have a lifetime history of fitness and who used to do more than 12 hours a week but do less now.
- Do not step off until "0000" disappears.

Measuring Weight Only

Press the $\triangle$ button. After 2 or 3 seconds, "0.00kg" (0.0lb, 0.0stlb) is displayed.

When "0.00kg" (0.0lb, 0.0stlb) is displayed, step onto the monitor.

The body weight measurement is displayed, and the monitor power turns off automatically after approximately 20 seconds.

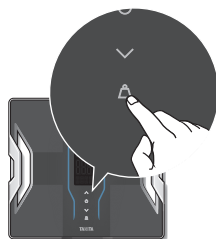

Body Composition Guide

Who Can Use a Body Composition Monitor?

This Body Composition Monitor is intended for adults aged 18-99 years. Children aged 5-17 years can use the monitor for weight, body fat percentage, muscle mass and bone mass: the other features are not applicable to children.

Thanks to major investment in the latest BIA Technology and sports science research, TANITA has upgraded the Athlete Mode function to make measurements more accurate and suitable for a wider range of users.

Who should use Athlete Mode?

Adults aged 18 years and over who either:

- Train or exercise for 12 hours or more a week and have been doing so for at least six months.
- Are body builders.
- Are professional athletes who want to monitor their progress at home.
- Have a lifetime history of fitness and who used to do more than 12 hours a week but do less now.

Pregnant women should only use the weight function.

Other functions are not intended for use when pregnant.

This Body Composition Monitor is intended for home use only. It is not intended for professional use in places such as hospitals or medical or fitness facilities. It is not designed for such heavy usage. Using the monitor in this type of professional environment will invalidate the warranty.

This Body Composition Monitor product provides readings for informational purposes only.


This product is not intended to diagnose or treat any disease or abnormalities. Please consult your physician if you have any questions or concerns related to your health.

How Does a Body Composition Monitor Work?

Tanita Body Composition Monitors calculate your body composition using Dual Frequency Bioelectrical Impedance Analysis (BIA). Safe, low-level electrical signals are passed through the body via the Tanita foot pads on the monitor platform. The signal can flow easily through fluids in muscles and other body tissue but meets resistance as it passes through body fat, because body fat only contains a small amount of fluid. This resistance is called impedance. The impedance readings are then entered into medically researched mathematical formulas to calculate your body composition.

The Tanita RD-901 incorporates medical grade Dual Frequency BIA technology, allowing you to have the highest body composition accuracy in the comfort of your home.

Research has shown that using two different bioelectrical impedance frequencies, provides essential data of a person's intracellular and extracellular status. This advanced technology allows greater accuracy when calculating body composition measurements.

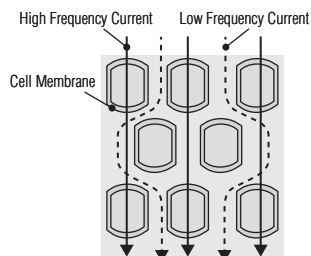

When Is the Best Time To Use My Body Composition Monitor?

Your body water levels naturally fluctuate throughout the day and night. Any significant changes in body water may affect your body composition readings. For example, the body tends to be dehydrated after a long night sleep so if you take a reading first thing in the morning your weight will be lower and your body fat percentage higher. Eating large meals, drinking alcohol, menstruation, illness, exercising, and bathing may also cause variations in your hydration levels.

To get the most reliable reading it is important to use your Body Composition Monitor at a consistent time of day under consistent conditions. A good time to take measurements is before your evening meal.

What Is Body Fat Percentage? (Applicable age 5-99)


Body fat percentage is the amount of body fat as a proportion of your body weight.

Reducing excess levels of body fat has been clinically proven to reduce the risk of certain conditions such as high blood pressure, heart disease, diabetes and cancer.

The chart below shows the healthy ranges of body fat.

Body Fat Ranges for Standard Children 1

Body Fat Ranges for Standard Adults 2,3,4

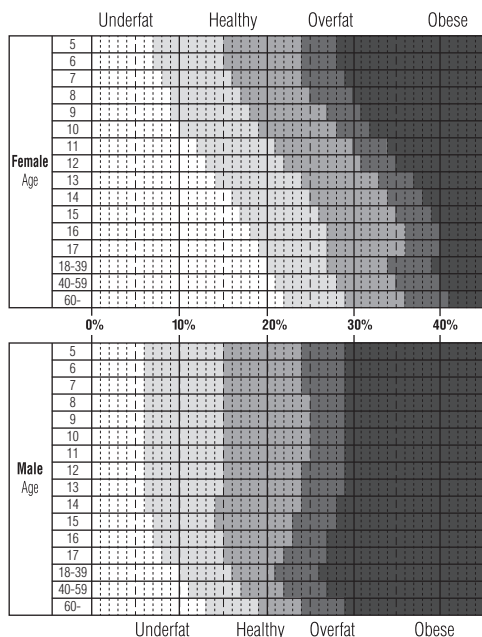

1. Body Fat Reference Curves for children Targeted at BMJ (British Medical Journal) Draft 1-AMP 19 June 2004 (by Dr Andrew)
 2. Gallagher D et al. Am J Clin Nutr 2000, 72:694-701 Healthy percentage body fat ranges: an approach for developing guidelines based on body mass index.
 3. Based on NIH/WHO BMI Guidelines.
 4. As reported by Gallagher, et al, at NY Obesity Research Center.
- To determine the percentage of body fat that is appropriate for your body, consult your physician.

Underfat: Below the healthy body fat range. Increased risk of health problems.

-Healthy • + Healthy: Within the healthy body fat percentage range for your age/gender.

Overfat: Above the healthy range. Increased risk of health problems.

Obese: Far above the healthy body fat range.

Greatly increased risk of obesity-related health problems.

What Is Total Body Water Percentage? (Applicable age 18-99)

Total Body Water Percentage is the total amount of fluid in a person's body expressed as a percentage of their total weight. Water plays a vital role in many of the body's processes and is found in every cell, tissue and organ. Maintaining a healthy total body water percentage ensures the body functions efficiently and reduces the risk of developing associated health problems.

Your body water levels naturally fluctuate throughout the day and night. Your body tends to be dehydrated after a long night's sleep and there are differences in fluid distribution between day and night. Eating large meals, drinking alcohol, menstruation, illness, exercising and bathing may cause variations in your hydration levels.

Your body water percentage reading should be used as a guide and should not be used to specifically determine your recommended total body water percentage. It is important to look for long-term changes in total body water percentage and maintain a consistent, healthy total body water percentage.

Drinking a large quantity of water in one sitting will not instantly change your body water level. In fact, it will increase your body fat reading due to the additional weight gain. Please monitor all readings over time to track relative changes.

Every individual is different, but as a guide the average total body water percentage ranges for a healthy adult are:

Female: 45 to 60%

Male: 50 to 65%

Source: Tanita Research

Note: The total body water percentage will tend to decrease as the percentage of body fat increases. A person with a high percentage of body fat may fall below the average body water percentage. As you lose body fat, the total body water percentage should gradually move towards the typical range given above.

What Is Visceral Fat Rating? (Applicable age 18-99)

This function indicates the amount of visceral fat in your body.

Visceral fat is the fat that is in the internal abdominal cavity, surrounding the vital organs in the abdominal area. Research shows that even if your weight and body fat remains constant, as you get older the distribution of fat changes and is more likely to shift to the abdominal area especially post menopause. Ensuring you have healthy levels of visceral fat may reduce the risk of certain diseases such as heart disease, high blood pressure, and the onset of type 2 diabetes.

The Tanita Body Composition Monitor provides you with a visceral fat rating from 1 – 59.

Rating 1 – 9.5

Indicates you have a healthy amount of visceral fat. Continue monitoring your rating to ensure that it stays within this healthy range.

Rating 10 – 59

Indicates you have an excess amount of visceral fat. Consider making changes in your lifestyle by changing your diet or exercising more.

Source : Data from Columbia University (New York) & Tanita Institute (Tokyo)

Note: You may have a high visceral fat level even if you have a low body fat rate.

- Consult a physician for an accurate medical diagnosis.

What Is Basal Metabolic Rate (BMR)? (Applicable age 18-99)

Your Basal Metabolic Rate (BMR) is the minimum level of energy your body needs when at rest to function effectively. This includes the functioning of your respiratory and circulatory organs, neural system, liver, kidneys, and other organs. You continue to burn calories even when sleeping.

About 70% of the calories your body uses every day are used by your basal metabolism. In addition, energy is used when doing any kind of activity. The more vigorous the activity, the more calories are burned. This is because skeletal muscle (which accounts for approximately 40% of your body weight) acts as your metabolic engine and uses a large amount of energy. Your basal metabolism is greatly affected by the amount of muscle you have, so increasing your muscle mass helps increase your basal metabolism.

By studying healthy individuals, scientists have found that people's metabolic rate change as they age. Basal metabolism rises as a child grows. After peaking at the age of around 16 or 17, it typically starts to gradually decrease.

Having a higher basal metabolism increases the number of calories used and helps decrease the amount of body fat. A low basal metabolic rate makes it harder to lose body fat and overall weight.

How Does a Tanita Body Composition Monitor Calculate BMR?

The basic method of calculating Basal Metabolic Rate (BMR) is a standard equation that uses weight and age.

Tanita has conducted in-depth clinical research into the relationship of BMR and body composition resulting in a much more accurate and personalized reading for the user based on impedance measurements. This method has been medically validated using indirect calorimetry (measuring breathcomposition).*

*Reliability of the equation for Basal Metabolic Rate: At 2002 Nutrition Week: Scientific and Clinical Forum and Exposition

Title: International Comparison: Resting Energy Expenditure Prediction Models:

The American Journal of Clinical Nutrition.

What Is Metabolic Age? (Applicable age 18-99)

This function calculates your BMR and indicates the average age associated with the type of metabolism.

If your metabolic age is higher than your actual age, it is an indication that you need to improve your metabolic rate. Increased exercise builds healthy muscle tissue, which improves your metabolic age.

The reading is shown as a number between 12 and 90. Values under 12 are displayed as "12" and over 90 are displayed as "90".

What Is Muscle Mass? (Applicable age 5-99)

This function indicates the weight of muscle in your body. The muscle mass displayed includes skeletal muscle, smooth muscle (such as cardiac and digestive muscle) and the water contained in these muscles.

Muscles play an important role as they act as an engine that consumes energy. As your muscle mass increases, your energy consumption increases helping you reduce excess body fat levels and lose weight in a healthy way.


What Is Physique Rating? (Applicable age 18-99)

This feature assesses your physique according to the ratio of body fat and muscle mass in your body.

As you become more active and reduce the amount of body fat, your physique rating also changes accordingly. Even if your weight does not change, your muscle mass and body fat levels may be changing to make you healthier and reduce your risk of certain diseases. Each person should set their own goal of which physique they want to achieve, and follow a diet and fitness program to meet that goal.

Result	Physique Rating	Explanation of Physique Rating Results
1	Hidden Obese	Small Frame Obese
		This person seems to have a healthy body type based on physical appearance. However, they have a high body fat percentage with low muscle mass level.
2	Obese	Medium Frame Obese
		This person seems to have a healthy body type based on physical appearance. However, they have a high body fat percentage with low muscle mass level.
3	Solidly-built	Large Frame Obese
		This person has both a high body fat percentage and a high muscle mass.
4	Under exercised	Low Muscle and Average Body Fat Percentage
		This person has an average body fat percentage and a less than average muscle mass level.
5	Standard	Ave. Muscle & Ave. Body Fat Percentage
		This person has average levels of both body fat and muscle mass.
6	Standard Muscular	High Muscle & Ave. Body Fat Percentage (Athlete)
		This person has an average body fat percentage and higher muscle mass level than average.
7	Thin	Low Muscle & Low Fat
		Both body fat percentage and muscle mass are lower than average.
8	Thin and muscular	Thin and muscular (Athlete)
		This person has a lower than average body fat percentage while having adequate muscle mass.
9	Very Muscular	Very Muscular (Athlete)
		This person has a lower than average body fat percentage while having above-average muscle mass.

Source: Data from Columbia University (New York) & Tanita Institute (Tokyo)

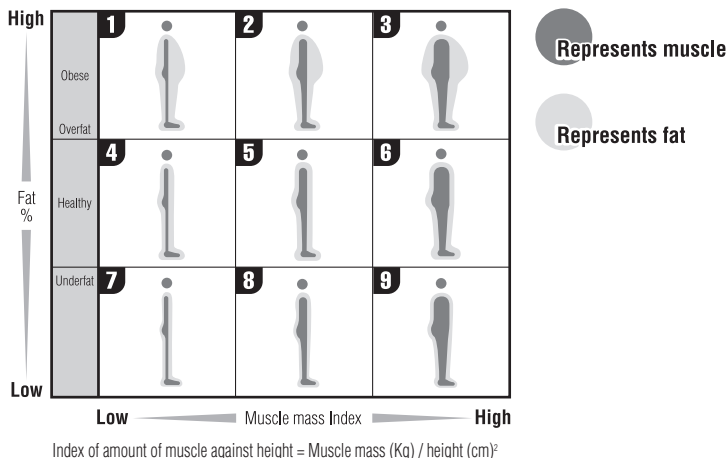

What Is Bone Mass? (Applicable age 5-99)

This function indicates the amount of bone (bone mineral level, calcium or other minerals) in the body.

Research has shown that exercise and the development of muscle tissue are related to stronger, healthier bones. While it is unlikely that there will be noticeable changes in bone structure over short periods, it is important that you develop and maintain healthy bones by having a balanced diet and plenty of exercise. People worried about bone disease should consult their physician. People who suffer from osteoporosis or low bone density due to age (young or old), pregnancy, hormonal treatment or other causes may not get accurate bone mass readings.

Below is the result of estimated bone masses of persons aged 20 to 40, who are said to have the largest amounts of bone mass by weight. (Source: Tanita Body Weight Science Institute)

Please use the charts below as a guide for comparing your bone mass readings.

Women: Average estimated bone mass

Men: Average estimated bone mass

Note:

- The following types of people may obtain varying results and should use the results for reference only.
 - Elderly people
 - Women during or after menopause
 - People undergoing hormone therapy
 - "Estimated bone mass" is a value estimated statistically based on its correlation with fat-free tissue (body tissue other than fat).
- "Estimated bone mass" does not give a direct judgment on the hardness or strength of bones or the risk of bone fractures. If you are concerned about your bones, it is recommended that you consult a specialist.

What Is Body Mass Index (BMI)? (Applicable age 5-99)

This function indicates the relationship between height and weight.

The World Health Organization (WHO) considers an index of 18.5 to 25 as optimal.

Troubleshooting

If the following problems occur, follow the instructions below.

“Lo” appears on the display.

Check the orientation of the batteries, and insert them correctly.

The batteries are running low. Replace the batteries promptly.

Refer to page 6.

“OL” appears on the display during use

The maximum weight (weighing capacity) 200kg is exceeded. (Measurements cannot be taken)

“Err” appears on the display after measurement

Are the personal data settings correct?

Did you select the wrong personal data number?

Measurement items have exceeded the measurement range. (Measurements cannot be taken)

Cannot measure body composition even though personal data is set

The personal data settings are not complete. You need to measure body composition once immediately after entering the personal data.

Cannot connect my iPhone to the Body Composition Monitor

“Err pair” or “Err DATA” appears on the display

Are batteries inserted correctly in the Body Composition Monitor?

Are the batteries running low?

Are the iPhone and body composition meter within communication range?

The communication range for The Body Composition Monitor is 5m in an unobstructed location.

Is the “Bluetooth” setting in the iPhone “Settings” turned on?

“Err pair” appears repeatedly.

Delete the pairing information from “Settings” > “Bluetooth” on the iPhone, and pair the iPhone with the Body Composition Monitor again.

“Err UUID” appears repeatedly.

Delete the pairing information from “Settings” > “Bluetooth” on the iPhone, and pair the iPhone with the Body Composition Monitor again.

“Err c9” appears on the display

Remove and replace the batteries from the body composition meter, and set the date and time again.

The date and time settings have been changed

When using the Body Composition Monitor linked with an iPhone, the date and time settings on the Body Composition Monitor are overwritten with those from the iPhone

The “Input Data” button in the app does not work

Is the “Bluetooth” setting in the iPhone “Settings” turned on?

Linking between the app and the Body Composition Monitor is only enabled when the Body Composition Monitor is turned off.

Check that the power of the Body Composition Monitor is turned off and try again.

Specifications

RD-901

Weight Capacity

200 kg (440 lb)(31st 6 lb)

Weight Increments

0-100 kg/0.05 kg 100-200 kg/0.1 kg
0-200 lb/0.1 lb 200-440 lb/0.2 lb 0-31st 6 lb/1 lb

Body Fat %

5-99 years old

Body Mass Index

5-99 years old

Muscle Mass

5-99 years old

Physique Rating

18-99 years old

Bone Mass

5-99 years old

Visceral Fat

18-99 years old

Basal Metabolic Rate

18-99 years old

Metabolic Age

18-99 years old

Body Water %

18-99 years old

Personal Data

4 memories

Power Supply

4×AA Alkaline batteries(included)

Power Consumption

250mA maximum

Measuring Current

50kHz, 6.25 kHz, 100μA

Communication Method

Bluetooth version 4.0 (Low Energy support)

Disposal



This is an electronic device.

Please dispose of it as an electronic device, not as general household waste. Please follow the regulations in your local region when disposing of this device.



⚠ Not allowed to mix batteries with consumer wastes!

As consumer you are legally bound to return used or discharged batteries. You can deposit your old batteries at the public collecting points in your town, or wherever the corresponding batteries are sold and specifically marked collecting boxes have been set up. In case of scrapping the apparatus, the batteries should be removed from it and deposited at the collecting points as well.



This product meets the following requirements;

1. EMC Directive
2. R&TTE Directive
3. RoHS Directive

•The product design and specifications may be changed at any time without prior notice.

•The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such mark by Tanita Corporation is under license.

•Apple, the Apple logo, and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

Manufacturer

TANITA Corporation

1-14-2, Maeno-Cho, Itabashi-ku,
Tokyo, Japan 174-8630
Tel: +81(0)-3-3968-7048
Fax: +81(0)-3-3968-2661
www.tanita.co.jp
ISO 9001 Certified

TANITA Corporation of America Inc.

2625 South Clearbrook Drive
Arlington Heights, IL 60005 USA.
Tel: +1-847-640-9241
Fax: +1-847-640-9261
www.tanita.com

Representative

TANITA Europe B.V.

Hoogoorddreef 56-E
1101 BE Amsterdam
The Netherlands
Tel: +31-20-560-2970
Fax: +31-20-560-2988
www.tanita.eu

TANITA Health Equipment H.K. Ltd.

Unit 301-303 3/F Wing On Plaza,
62 Mody Road, Tsimshatsui East,
Kowloon, Hong Kong
Tel: +852-2834-3917
Fax: +852-2838-8667
www.tanita.asia

TANITA (Shanghai) Trading Co. Ltd.

Room 8005, 887 Huai Hai Zhong Lu, Shanghai
The People's Republic of China
Tel: +86-21-6474-6803
Fax: +86-21-6474-7901
www.tanita.com.cn

TANITA India Private Limited

A-502, Mittal Commercial,
Off. M. V. Road (Andheri Kurla Road),
Marol, Andheri - East,
Mumbai, India 400 059
Tel: +91-22-3192-6107
Fax: +91-22-2859-9143
www.tanita.co.in