

TANITA

Consumer Product Catalog



WORLD LEADING BIA Brand

www.tanita.asia

TANITA

Trusted by Professionals for Over 70 Years

Founded in 1944, TANITA is the worldwide leading brand of medical, commercial and consumer body fat monitors/scales. In 1992, TANITA introduced the World's First body fat scale for personal use. Over the years, TANITA's accuracy, innovation and excellent product quality have set a golden standard for Bio-electrical Impedance Analysis (BIA) monitors and are well trusted by physicians, clinicians and health & fitness experts worldwide. TANITA aims to develop unique health management solution in all-rounded aspects that enable people to enjoy fuller and healthier lives.



TANITA Headquarter in Tokyo, Japan



TANITA Café in Tokyo, Japan



TANITA FITME in Tokyo, Japan

All-rounded Healthy Habits for Happiness



Professional BIA Technology

- Base on Bio-Electrical Impedance Analysis (BIA) to estimate the amount of body fat, muscle and other readings
- Send electrical signal throughout the body via metal footplates (and through handgrips for segmental monitors)
- Passes freely through blood, water and fluid within muscle tissue, but meets resistance when contacting with body fat
- Accurately measure body composition readings with reference to a scientific equation reference to users' gender, height and other information

Accurate and Quality Assurance

- At least four footplates are required for body composition monitor to ensure maximum accuracy of body composition readings
- Incorporate sound independent medical research spanning 25 years with 100,000 people worldwide in its studies, to ensure maximum accuracy of readings
- Estimate results based on electrical impedance at 50kHz frequency, which is the most reliable frequency referenced to DEXA
- Consumer products are certified by Japan JIS standard

TANITA HealthPlanet app

- Instantly record your results and track your progress.
- Easy-to-read charts and graphs for monitoring physical changes and body changes across different period of time.
- Set your own targets and achieve your peak performance.



Visualize your body composition for better motivation and insights




Overview



Segmental Muscle / Body Fat Readings (available for RD-545)



Detailed list of body composition reading

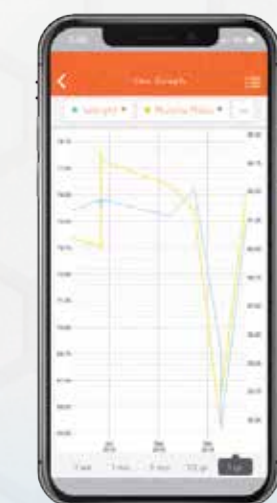

Graphs showing the body composition trend

TANITA Body Composition Measurements



Body Fat Ranges (Adults and Children)

Body Fat % is the amount of fat as a proportion of body weight. Body fat is vital to basic bodily functions such as regulating body temperature, storing vitamins and cushioning joints. Yet, too much fat can damage your health. The indicator shows if the amount of body fat falls in the Underfat, Healthy, Overfat or Obese categories.



Total Body Water %

Total Body Water % is the total amount of fluid in the body expressed as a percentage of total body weight. Being well hydrated is important to support you to have good sports performance and well-being.



Muscle Quality Score

Muscle Quality Score indicates the condition (quality) of muscle which changes according to factors like age and exercise level. The score is rated at a full score of 100 and classified into 3 categories: High, Average and Low.



Daily Calories Intake (DCI)

Daily Calories Intake automatically calculates the number of calories an individual should consume within 24 hours in order to maintain their current weight.



Basal Metabolic Rate (BMR)

Basal Metabolic Rate (BMR) is the daily minimum level of energy or calories you need at rest for your body to function effectively. A person at rest with a high BMR can burn more calories than person with a low BMR.



Metabolic Age

Metabolic Age indicates what age level the user's metabolism is associated with.



Muscle Mass

The weight of muscle in your body. Muscle plays an important role because they act as the body's engine in consuming energy (calories). As you exercise more, your muscle builds up, which in turn means increasing ability to burn more calories.



Visceral Fat Rating

This rating indicates the level of internal fat in the abdomen that is scientifically proved to increase the risk of diabetes, hypertension and other metabolic syndromes. It is shown in a rating of 1 to 59.



Physique Rating

Physique Rating is an assessment of your physique according to the ratio of body fat and muscle mass in your body. With increasing activity level, weight may not change but the ratio of body fat and muscle may alter, which will change the user's overall physique.

Physique Rating			
Fat	Obese	Over Fat	Healthy
	Hidden obese	Obese	Solidly-Build
	Under Exercised	Standard	Standard Muscular
	Under Fat	Thin	Thin and Muscular
	Thin	Thin and Muscular	Very Muscular
	-	0	+
	Muscle Mass		

Physique Rating will be displayed in a range of 1-9 as showing in the following graph:

Index of amount of
muscle against height =
muscles mass (kg) / height (cm²)



Athlete Mode

TANITA body composition monitors incorporate the Athlete Mode for anyone who has carried out more than 12 hours of intensive aerobic exercise a week for at least 6 months and has a resting heartbeat of under 60bpm. It also applies to people who have a lifetime history of fitness - who used to do more than 12 hours exercise a week but do less now, still classified to athlete type.



Bone Mass

Bone Mass is the weight of bone mineral (including calcium or other minerals) in your body. It is useful to monitor and check for significant changes over time.

*TANITA's body composition monitors do not measure bone density. Please consult a physician if you are concerned about osteoporosis.



Segmental Measurements

With the 8-electrode segmental technology, monitors show individual's body composition measurements for each arm, leg and trunk area. Essential for the dedicated sports person to assess the impact of their training routine on their bodies.



FitPLUS feature for Women

After user keys in the start day of her menstrual cycle, the feature will hint her about different phases of the cycle to help maximize calories burnt and minimize fatigue.

Red light : Relaxation days

Amber light : Ovulation

Green light : Optimal active days



MBA Judgment


The MBA (My Body Analyzer) judgment is an evaluation based on your muscle condition, muscle fibers properties and muscle balance. This score benchmark your body ability by comparing with our 1,000+ professional athlete database.




Dual Technology

Professional Dual Frequency

- TANITA RD-series features the pioneering medical grade **Dual Frequency BIA Technology**
- Using two different bio electrical impedance frequencies provides essential data of a person's intra cellular and extra cellular status
- Provides comprehensive and precise body composition results as well as one's muscle condition – **"Muscle Quality"**

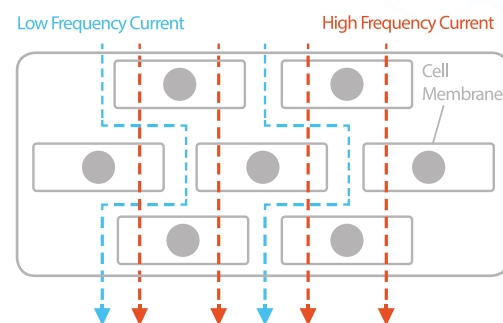


Muscle Quality Score

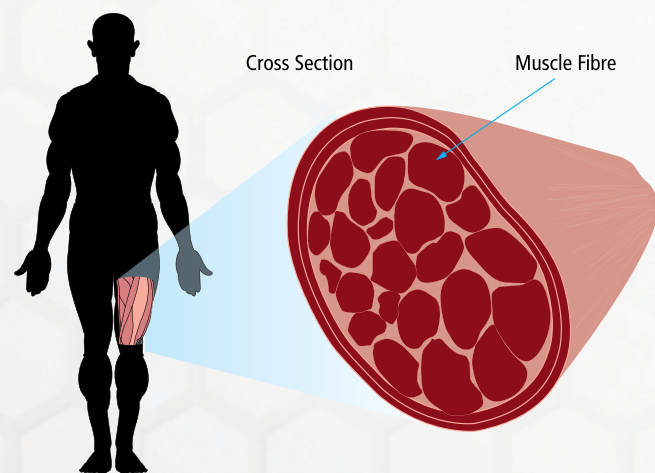

Muscle Strength is mainly related to Muscle Mass and Muscle Quality. Both aging and exercise level would change the muscle condition, i.e. the Muscle Quality.

Monitoring this Score can give you valuable insights on how to customize your training program to boost performance and avoid injuries.

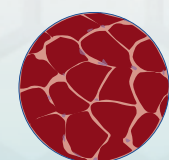


Muscle Quality is **HIGHER**

Higher density of muscle fibre with little fat, water and connective tissue, etc.



Maintain regular exercise to ensure your Muscle Quality stays within the healthy

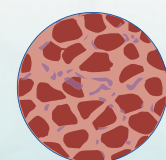

Muscle Quality is **LOWER**

Muscle fibre is thinner and fewer as compared with other substances.



Increase exercise on a regular basis to improve your Muscle Quality.

RD-545



Dual Technology



Operation Video



Dual Frequency Segmental Body Composition Monitor for Smartphones

The First & the Best

- Accurate – Pioneering Dual Frequency BIA technology for professional accuracy.
- In-depth – Thorough 26 whole-body/ segmental analysis obtained in less than 60 seconds.
- Unique in the market.
- Convenient – Access your body composition readings anytime, anywhere via TANITA HealthPlanet app.
- Made in Japan Quality Guarantee.


Technical Specifications

Max. Capacity/ Graduation	0 - 100 kg/ 0.05kg 100 - 200kg/ 0.1kg
Product Size (mm)	348 x 320 x 57mm
Product Weight	2.8kg
Batteries	AA x 4
No. of User memories	4
Display Language	English
Device Connection	Bluetooth
Supported OS	iOS 10 or later, Android 5.1 later

Packaging Information

Gift Box Size	364 x 353 x 75mm
Carton Size	381 x 318 x 368mm
Carton Quantity	4pcs/ carton
Carton Weight	14.3kg


Body Composition Measurements:

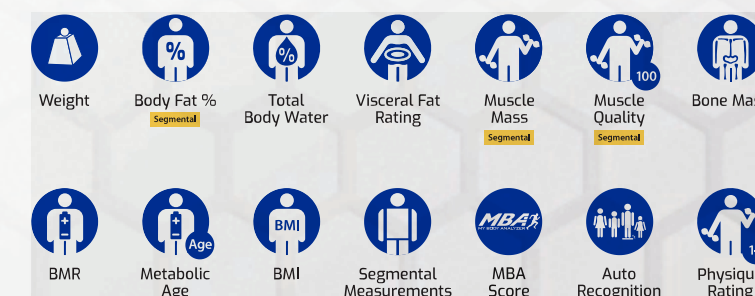

RD-953



Dual Technology

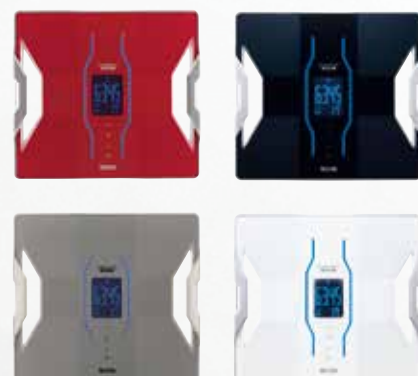

Operation Video



Dual Frequency Body Composition Monitor for Smartphones

- Dual frequency technology allows comprehensive body composition readings
- Quick data connection to free mobile app, allows you to review & track your body composition level anytime and anywhere
- 200kg max weight capacity with 50g graduation


Technical Specifications

Max. Capacity/ Graduation	0 - 100 kg / 0.05 kg 100 – 200 kg / 0.1 kg
Product Size (mm)	298 x 328 x 32mm
Product Weight	2.1kg
Batteries	AA x 4
No. of User memories	4
Display Language	English
Color Option	Black, White, Red, Grey
Device Connection	Bluetooth

Packaging Information

Gift Box Size	359 x 379 x 49mm
Carton Size	265 x 391 x 377mm
Carton Quantity	5pcs / carton
Carton Weight	14.3kg

Body Composition Measurements:

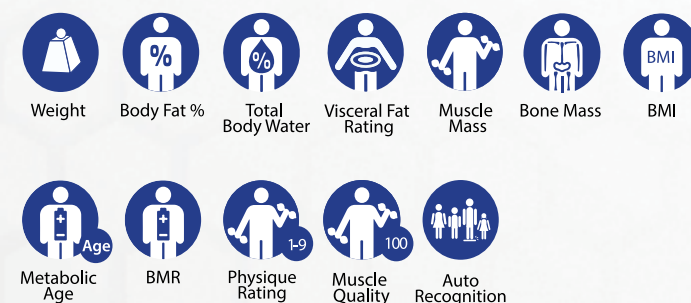

BC-402



Operation Video



Body Composition Monitor for Smartphones

- Apps to keep track your body changes
- Easy-to-clean glass platform and modern design to fit all homes
- One scale to store up to 5 family members data


Technical Specifications

Max. Capacity/ Graduation	150kg/ 0.1kg
Product Size (mm)	217 x 316 x 27mm
Product Weight	1.5kg
Batteries	AAA x 4
No. of User memories	5
Display Language	English
Color Option	Black, White
Device Connection	Bluetooth
Supported OS	iOS 10 or later, Android 5.1 later

Packaging Information

Gift Box Size	261 x 356 x 43mm
Carton Size	233 x 368 x 278mm
Carton Quantity	5pcs / carton
Carton Weight	9.5kg

Body Composition Measurements:

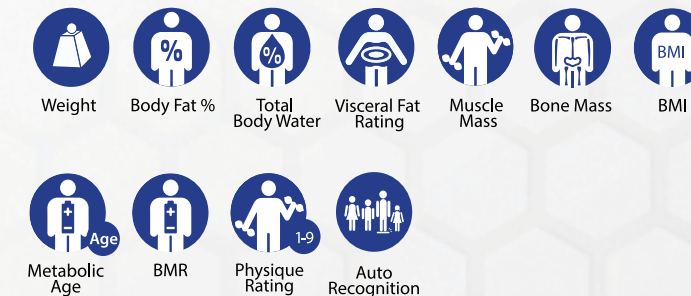

BC-333L




Body Composition Monitor for Smartphones

- Review and track your body changes via TANITA HealthPlanet app
- When measuring yourself, it also shows an user number that has not been measured for a certain period
- Extra-large display with showing battery power level
- 50g graduation
- One scale to store up to 5 user members data

Technical Specifications	
Max. Capacity/ Graduation	0 – 100kg / 0.05kg 100 – 150kg / 0.1kg
Product Size (mm)	274 X 310 X 36mm
Product Weight	1.5kg
Batteries	AA x 4
No. of User memories	5
Display Language	Japanese
Color Option	Pink, White
Device Connection	Bluetooth
Supported OS	iOS 10 or later, Android 5.1 later

Packaging Information	
Gift Box Size	326 x 323 x 44mm
Carton Size	330 x 240 x 338mm
Carton Quantity	5pcs / carton
Carton Weight	9.1kg

Body Composition Measurements:

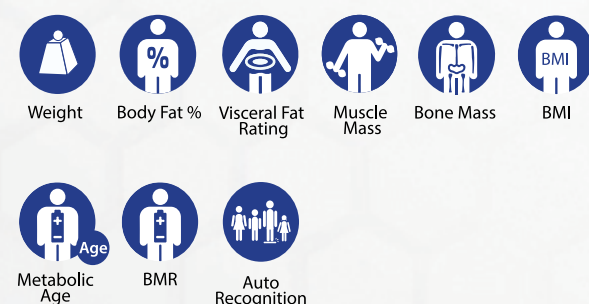

BC-545N



Operation Video



BC-313




Segmental Body Composition Monitor with Easy-To-Read Analytical Graphic Display

- High resolution color display including newly improved easy-to-read analytical graphics
- Measurement tracking on a weekly, fortnightly and annual basis, allows you to review your fitness progress instantly
- Segmental readings for fat and muscle displayed in 9 levels, easy to assess the distribution of different parts and corresponding weight in one glance

Technical Specifications		Packaging Information	
Max. Capacity / Graduation	150kg/0.1kg	Gift Box Size	378 x 454 x 69mm
Product Size (mm)	320 x 384 x 57mm	Carton Size	465 x 298 x 398mm
Product Weight	2.5kg	Carton Quantity	4pcs/ carton
Batteries	AA x 4	Carton Weight	14.4kg
No. of User memories	5		
Display Language	English		

Body Composition Measurements:

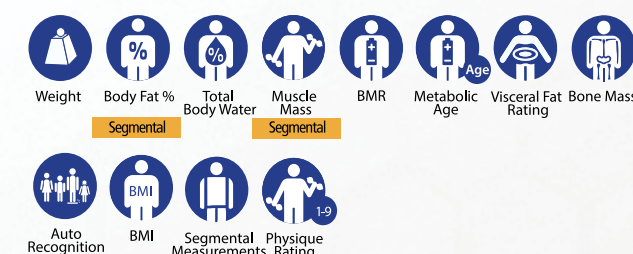


Glass Body Composition Monitor with Extra-Large Display and 50g Accuracy

- 50g graduation
- 29mm easy-to-read extra-large display
- Tempered glass surface with special plastic rim protection

Technical Specifications		Packaging Information	
Max. Capacity / Graduation	150kg/ 0.05kg	Gift Box Size	261 x 356 x 43mm
Product Size (mm)	217 x 316 x 27mm	Carton Size	233 x 368 x 278mm
Product Weight	1.5kg	Carton Quantity	5pcs/ carton
Batteries	AAA x 4	Carton Weight	9.5kg
Size of Display	29mm		
No. of User memories	5		
Display Language	English		
Color Option	Black, White, Red		

Body Composition Measurements:

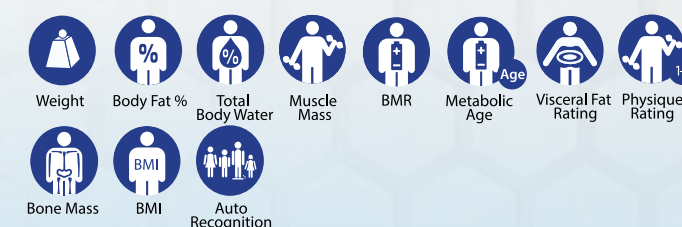

BC-G02




Lightweight Body Composition Monitor with FitPlus Feature

- Equipped with target weight setting and FitPlus feature, perfect companion for female
- The monitor will notify the user if there is a significant change in her body weight

Technical Specifications	
Max. Capacity / Graduation	150kg/ 0.1kg
Product Size (mm)	226 x 266 x 34mm
Product Weight	0.9kg
Batteries	AA x 4
No. of User memories	2
Size of Display	21.3mm
Display Language	English
Color Option	White, Lime Green, Pink

Packaging Information	
Gift Box Size	292 x 300 x 56mm
Carton Size	270 x 330 x 320mm
Carton Quantity	4pcs/ carton
Carton Weight	8.3kg

Body Composition Measurements:

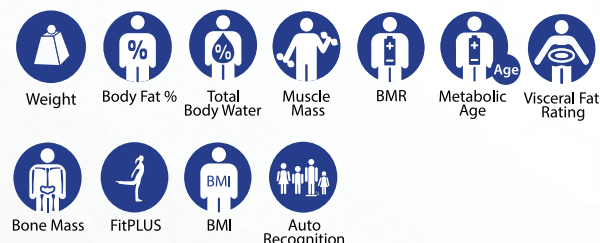

BC-541N




Lightweight Plastic Body Composition Monitor

- 30mm easy-to-read extra-large display
- Blue Backlit LCD Display

Technical Specifications	
Max. Capacity / Graduation	0 – 150kg/ 0.1 kg
Product Size (mm)	285 x 278 x 36 mm
Product Weight	1.2kg
Batteries	AA x 4
No. of User memories	5
Size of Display	30mm
Display Language	English
Color Option	Black, White

Packaging Information	
Gift Box Size	334 x 361 x 52mm
Carton Size	240 x 370 x 365mm
Carton Quantity	4pcs/ carton
Carton Weight	7kg

Body Composition Measurements:

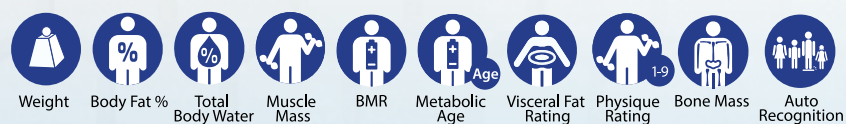

Operation Video



BC-730



Compact and Lightweight Body Composition Monitor

- 880g lightweight and compact

Technical Specifications	
Max. Capacity / Graduation	150kg/ 0.1kg
Product Size (mm)	216 x 260 x 33mm
Product Weight	880g
Batteries	AA x 4
No. of User memories	5
Size of Display	21.3mm
Display Language	English
Color Options	Pink, Green, White

Packaging Information	
Gift Box Size	253 x 300 x 44mm
Carton Size	316 x 288 x 274mm
Carton Quantity	6pcs/ carton
Carton Weight	7.2kg

Body Composition Measurements:



BC-731




Compact Glass Body Composition Monitor

- A4 compact size with large, easy-to-read display
- Tempered glass surface, easy to clean

Technical Specifications	
Max. Capacity / Graduation	150kg/ 0.1kg
Product Size (mm)	210 x 297 x 25mm
Product Weight	1.4kg
Batteries	AAA x 4
No. of User memories	5
Display Language	English
Color Option	Pink, White, Black

Packaging Information	
Gift Box Size	254 x 337 x 41mm
Carton Size	223 x 349 x 271mm
Carton Quantity	5pcs/ carton
Carton Weight	8.3kg

Body Composition Measurements:

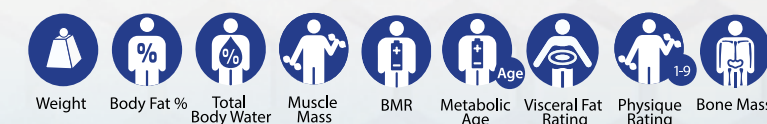

Operation Video




Compact Glass Body Composition Monitor

- “Active Level” indicates leg muscle level
- Easy-to-read arrow indicator showing you the trend compared with previous result

Technical Specifications	
Max. Capacity / Graduation	0 – 100kg/ 0.1kg, 100 – 150kg/ 0.2kg
Product Size (mm)	215 x 259 x 33mm
Product Weight	1kg
Batteries	AAA x 4
No. of User memories	5
Display Language	Japanese
Color Option	Red, White

Packaging Information	
Gift Box Size	293 x 337 x 45mm
Carton Size	245 x 350 x 311mm
Carton Quantity	5pcs/ carton
Carton Weight	7.1kg

Body Composition Measurements:

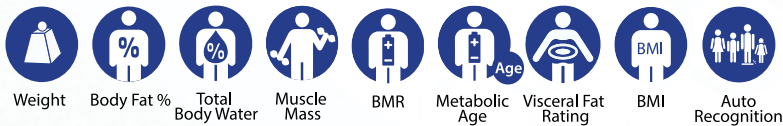

BC-543




Glass Body Composition Monitor with Touch Buttons

- Tempered Glass surface with easy-to-operate touch buttons

Technical Specifications	
Max. Capacity / Graduation	150kg/ 0.1kg
Product Size (mm)	320 x 300 x 44mm
Product Weight	2.4kg
Batteries	AA x 4
No. of User memories	4
Size of Display	32mm
Display Language	English

Packaging Information	
Gift Box Size	370 x 365 x 70mm
Carton Size	378 x 293 x 385mm
Carton Quantity	4pcs/ carton
Carton Weight	14.3kg

Body Composition Measurements:

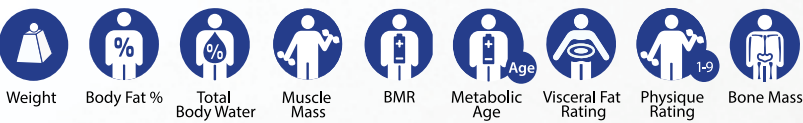

High Capacity Glass Body Composition Monitor

- Features high weight capacity and extra wide platform
- Results are clearly displayed on an extra-large LCD screen with backlight

Technical Specifications	
Max. Capacity / Graduation	200kg/ 0.1kg
Product Size (mm)	300 x 340 x 43mm
Product Weight	2.7kg
Batteries	AA x 4
No. of User memories	4
Size of Display	32mm
Display Language	English

Packaging Information	
Gift Box Size	374 x 390 x 68mm
Carton Size	306 x 410 x 390mm
Carton Quantity	4pcs/ carton
Carton Weight	14kg

Body Composition Measurements:



BC-541




Plastic Body Composition Monitor with Touch Buttons

- Durable grey plastic platform
- Easy-to-operate touch buttons

Technical Specifications	
Max. Capacity / Graduation	150kg / 0.1kg
Product Size (mm)	302 x 307 x 43.2mm
Product Weight	2.4kg
Batteries	AA x 4
No. of User memories	4
Size of Display	32mm
Display Language	English

Packaging Information	
Gift Box Size	355 x 355 x 64mm
Carton Size	365 x 273 x 366mm
Carton Quantity	4pcs/ carton
Carton Weight	12.6kg

Body Composition Measurements:

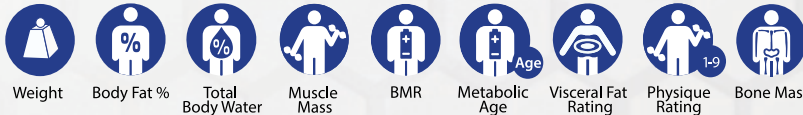

BC-758



Compact Glass Body Composition Monitor

- Tempered glass surface with plastic rim
- 29mm large text display
- Easy-to-step-on wide weighing platform
- Easy-to-read arrow indicator showing the trend compared with previous readings

Technical Specifications	
Max. Capacity / Graduation	0 – 100kg/ 0.1kg, 100 – 150kg/ 0.2kg
Product Size (mm)	293 x 286 x 25mm
Product Weight	1.7kg
Batteries	AAA x 4
No. of User memories	5
Display Language	Japanese
Color Option	Night Blue, Rose Pink, White

Packaging Information	
Gift Box Size	335 x 344 x 41mm
Carton Size	228 x 347 x 362mm
Carton Quantity	5pcs/ carton
Carton Weight	11.1kg

Body Composition Measurements:



BC-760




Plastic Body Composition Monitor

- Arrow indicator available for comparison with last measurement

Technical Specifications	
Max. Capacity / Graduation	150kg / 0.1kg
Product Size (mm)	215 x 259 x 33mm
Product Weight	1kg
Batteries	AA x 4
No. of User memories	5
Size of Display	27mm
Display Language	Japanese
Color Option	Red, Brown, White

Packaging Information	
Gift Box Size	293 x 337 x 45mm
Carton Size	245 x 350 x 311mm
Carton Quantity	5pcs / carton
Carton Weight	7.1kg

Body Composition Measurements:

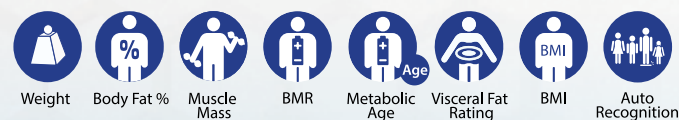

BC-759



Operation Video



BC-G01




Compact & Lightweight Body Composition Monitor

- Lightweight - 880g lightweight, easy to carry around
- Compact – Smaller than A4 size

Technical Specifications	
Max. Capacity / Graduation	0 – 100kg/ 0.1kg, 100 – 150kg/ 0.2kg
Product Size (mm)	216 x 260 x 33mm
Product Weight	880g
Batteries	AAA x 4
No. of User memories	5
Display Language	Japanese
Color Option	White, Light Blue, Pink, Brown

Packaging Information	
Gift Box Size	253 x 300 x 44mm
Carton Size	316 x 288 x 274mm
Carton Quantity	6pcs/ carton
Carton Weight	7.2kg

Body Composition Measurements:

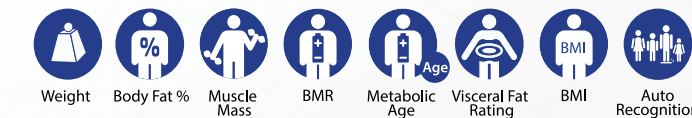


Plastic Body Fat Monitor

Technical Specifications	
Max. Capacity / Graduation	150kg/ 0.1kg
Product Size (mm)	284 x 284 x 39mm
Product Weight	1.2kg
Batteries	AA x 4
No. of User memories	5
Size of Display	21.3mm
Display Language	English
Color Option	Dark Blue, White

Packaging Information	
Gift Box Size	295 x 304 x 47mm
Carton Size	309 x 317 x 312mm
Carton Quantity	6pcs/ carton
Carton Weight	8.7kg

Body Composition Measurements:

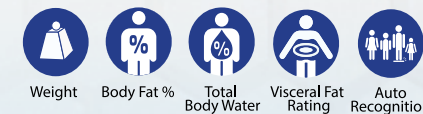

BC-770



Lightweight Plastic Body Composition Monitor

- 30mm easy-to-read extra-large display

Technical Specifications	
Max. Capacity / Graduation	0 - 100kg/ 0.1kg 100 - 150kg/ 0.2kg
Product Size (mm)	210 x 297 x 25mm
Product Weight	1.1kg
Batteries	AA x 4
Display Language	Japanese
Color Option	Dark Blue, Dark Grey

Packaging Information	
Gift Box Size	302 x 321 x 47mm
Carton Size	302 x 333 x 320mm
Carton Quantity	6pcs/ carton
Carton Weight	8.4kg

Body Composition Measurements:



NEW

HD-665



NEW

Carry Scale

- "Carry Mode" special calculator function to keep tracking baby's weight change. Let you easily review baby growth status through a graph
- Measure how much you feed for baby. Suitable for breast-feeding mothers and newborn baby
- 50g accuracy with automatic step-on activation
- Applicable for luggage weight measurement

Technical Specifications	
Max. Capacity/ Graduation	0 - 100 kg / 0.05kg 100 - 150kg / 0.1kg
Product Size (mm)	217 x 316 x 27mm
Product Weight	1.5kg
Batteries	AA x 4
Display Language	English
Color Option	White, Black

Packaging Information	
Gift Box Size	261 x 356 x 43mm
Carton Size	233 x 368 x 278mm
Carton Quantity	5pcs / carton
Carton Weight	9.5kg


Operation Guide



1. Press "Carry"




2. Measure the weight with the carried item

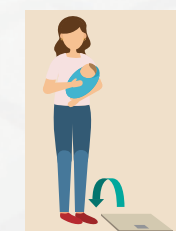

3. Step down




4. Measure the weight without the item

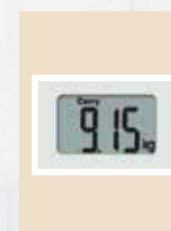

5. Calculation automatically shows the weight of the carried item

CA-100

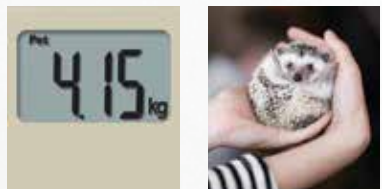

Pet Scale

- “Carry Mode” special calculator function to measure the weight of pet
- Setup regular training if your pet is overweight
- Feed the proper amount to avoid overfeeding
- 50g accuracy with automatic step-on activation


Technical Specifications

Max. Capacity/ Graduation	0 - 100 kg / 0.05kg 100 - 150kg / 0.1kg
Product Size (mm)	217 x 316 x 27mm
Product Weight	1.5kg
Batteries	AA x 4
Size of Display	36mm
Display Language	English

Packaging Information

Gift Box Size	261 x 356 x 43mm
Carton Size	233 x 368 x 278mm
Carton Quantity	5pcs / carton
Carton Weight	9.5kg

Operation Guide

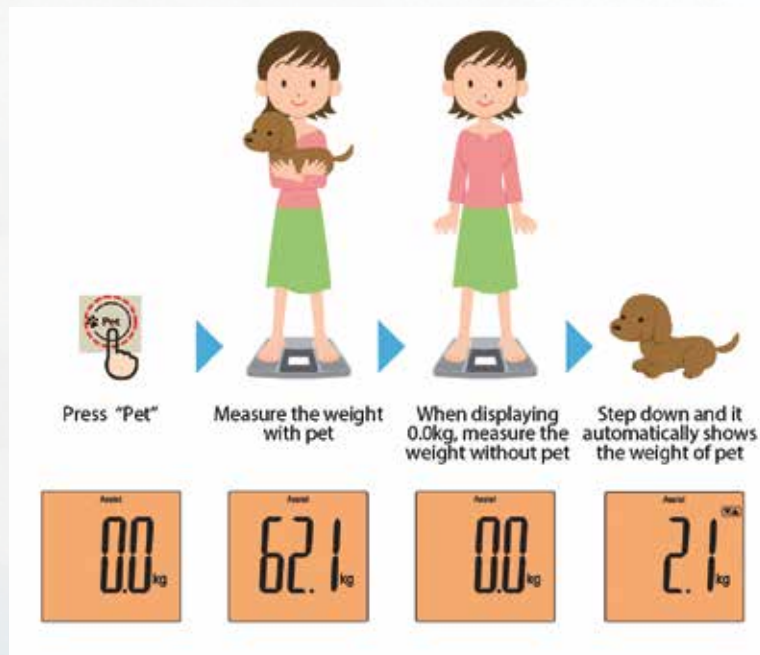

NEW

HD-366



HD-380



Digital Bathroom Scale with High Capacity and Blue Backlit Reverse LCD Display

Welcome by 5-star hotels

- Features extra-large (61mm) text display.

Technical Specifications

Max. Capacity/ Graduation	200kg / 0.1kg
Product Size (mm)	320 x 350 x 26mm
Product Weight	3kg
Batteries	AAA x 4
Size of Display	61mm

Packaging Information

Gift Box Size	421 x 428 x 50mm
Carton Size	203 x 433 x 440mm
Carton Quantity	4pcs/ carton
Carton Weight	15.3kg

Ultra Slim Glass Digital Bathroom Scale

Technical Specifications

Max. Capacity/ Graduation	150kg / 0.1kg
Product Size (mm)	300 x 320 x 22mm
Product Weight	1.7kg
Batteries	Lithium CR2032 x 1
Size of Display	30mm
Display Language	English
Color Option	Black, White, Blue, Light Pink, Light Blue, Pearl White, Pearl White, Black Square, Red, Orange, Green

Packaging Information

Gift Box Size	315 x 336 x 29mm
Carton Size	196 x 350 x 330mm
Carton Quantity	6pcs/ carton
Carton Weight	11.6kg

HD-381



Ultra Slim Digital Glass with Extra Large LCD Display

Technical Specifications	
Max. Capacity/ Graduation	150kg / 0.1kg
Product Size (mm)	300 x 320 x 22mm
Product Weight	1.7kg
Batteries	Lithium CR2032 x 1
Size of Display	52mm
Display Language	English
Color Option	Green, Red, Light Green

Packaging Information	
Gift Box Size	315 x 336 x 29mm
Carton Size	196 x 350 x 330mm
Carton Quantity	6pcs / carton
Carton Weight	11.6kg

HD-383




Digital Glass Scale with BMI Calculator

- After choosing your height using the dial button, the BMI reading and health indicator will come up immediately when stepping on the scale.

Technical Specifications	
Max. Capacity/ Graduation	150kg / 0.1kg
Product Size (mm)	300 x 320 x 33.7mm
Product Weight	2.2kg
Batteries	AAA x 2
Size of Display	52mm
Display Language	English
Color Option	Brown, Black, Green

Packaging Information	
Gift Box Size	324 x 336 x 45mm
Carton Size	292 x 350 x 343mm
Carton Quantity	6pcs / carton
Carton Weight	15.5kg

HD-382

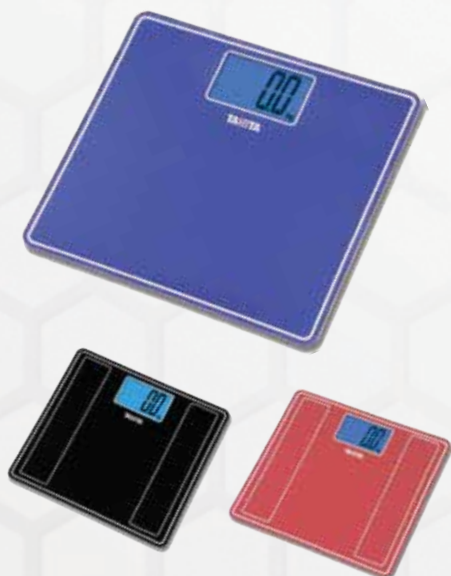


Digital Bathroom Scale with Extra Large

Technical Specifications	
Max. Capacity/ Graduation	150kg/ 0.1kg
Product Size (mm)	300 x 320 x 24mm
Product Weight	1.7kg
Batteries	AAA x 2
Size of Display	52mm
Display Language	English
Color Options	Metallic Black, Metallic Red, Metallic Blue

Packaging Information	
Gift Box Size	315 x 336 x 29mm
Carton Size	196 x 350 x 330mm
Carton Quantity	6pcs/ carton
Carton Weight	11.6kg

HD-384

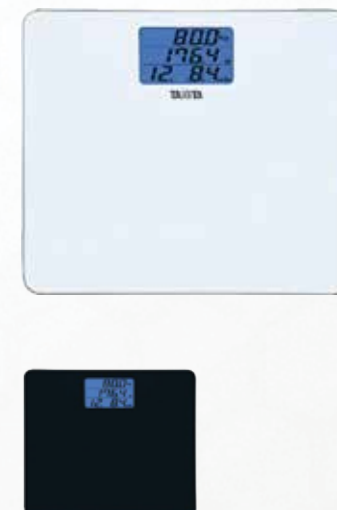

Digital Glass Scale with Simultaneous Multi Measurement Display

- Simultaneously displaying measurement in kg, lb and st-lb
- Features extra-large (52mm) blue LCD backlit display

Technical Specifications		Packaging Information	
Max. Capacity / Graduation	200kg / 0.1kg	Gift Box Size	421 x 428 x 50mm
Product Size (mm)	320 x 350 x 26mm	Carton Size	198 x 432 x 440mm
Product Weight	3kg	Carton Quantity	4pcs/ carton
Batteries	AAA x 4	Carton Weight	15.3kg
Size of Display	52mm		
Display Language	English		
Color Option	Black, White		



Digital Bathroom Scale with Extra Large Weighing Platform and High Capacity

Technical Specifications	
Max. Capacity / Graduation	200kg/ 0.1kg
Product Size (mm)	319 x 379 x 28mm
Product Weight	2.7kg
Batteries	AAA x 2
Size of Display	52mm
Display Language	English
Color Option	Black, Cream, Red, Pink, White

Packaging Information	
Gift Box Size	366 x 416 x 48mm
Carton Size	209 x 430 x 383mm
Carton Quantity	4pcs/ carton
Carton Weight	12.9kg



HD-395

Digital Bathroom Scale with Extra-Large Red LED Display

Fit for senior

- Features 54mm Extra Large Display

Technical Specifications	
Max. Capacity/ Graduation	150kg / 0.1kg
Product Size (mm)	278 x 279 x 34mm
Product Weight	1.1kg
Batteries	AA x 4
Size of Display	34mm
Display Language	English

Packaging Information	
Gift Box Size	314 x 319 x 45mm
Carton Size	290 x 331 x 332mm
Carton Quantity	6pcs/ carton
Carton Weight	8.4kg



Digital Bathroom Scale with Large Red LED Display

Technical Specifications	
Max. Capacity/ Graduation	150kg/ 0.1kg
Product Size (mm)	250 x 310 x 32mm
Product Weight	0.84kg
Batteries	AAA x 4
Size of Display	34mm
Display Language	English
Color Option	Wine Red, Black, White, Pink

Packaging Information	
Gift Box Size	286 x 351 x 40mm
Carton Size	260 x 363 x 305mm
Carton Quantity	6pcs/ carton
Carton Weight	7kg



HD-660

Compact Plastic Digital Bathroom Scale

Technical Specifications	
Max. Capacity/ Graduation	150kg/ 0.1kg
Product Size (mm)	200 x 244 x 33mm
Product Weight	0.7kg
Batteries	AA x 2
Size of Display	30mm
Display Language	English
Color Option	Pink, Purple, White, Blue

Packaging Information	
Gift Box Size	218 x 261 x 44mm
Carton Size	288 x 276 x 230mm
Carton Quantity	6pcs/ carton
Carton Weight	5.9kg



Plastic Digital Bathroom Scale

Technical Specifications	
Max. Capacity/ Graduation	150kg / 0.1kg
Product Size (mm)	280 x 280 x 35mm
Product Weight	0.99kg
Batteries	AA x 2
Size of Display	30mm
Display Language	English

Packaging Information	
Gift Box Size	300 x 328 x 43mm
Carton Size	279 x 343 x 318mm
Carton Quantity	6pcs/ carton
Carton Weight	8.2kg



Matt Finished Digital Bathroom Scale with Extra Large Blue Backlit LCD Display

Technical Specifications	
Max. Capacity/ Graduation	150kg/0.1kg
Product Size (mm)	288 x 288 x 35mm
Product Weight	1.1kg
Batteries	AA x 2
Size of Display	52mm
Display Language	English
Color Option	Black, White

Packaging Information	
Gift Box Size	300 x 328 x 43mm
Carton Size	279 x 343 x 318mm
Carton Quantity	6pcs / carton
Carton Weight	8.6kg



Solar Powered Digital Scale

Technical Specifications	
Max. Capacity/ Graduation	3 – 100kg/ 0.1kg, 100 – 150kg/ 0.2kg
Product Size (mm)	306 x 330 x 48mm
Product Weight	2.3kg
Batteries	NIL
Size of Display	24mm
Display Language	English

Packaging Information	
Gift Box Size	360 x 350 x 60mm
Carton Size	270 x 360 x 380mm
Carton Quantity	4pcs/ carton
Carton Weight	11.4kg



Mechanical Scale with Magnifying Window Display

Technical Specifications	
Max. Capacity/ Graduation	130kg/ 1kg
Product Size (mm)	252 x 254 x 43.5mm
Product Weight	1.4kg
Color Option	White, Black, Animal Girl, Tennis, Band

Packaging Information	
Gift Box Size	257 x 262 x 50mm
Carton Size	320 x 274 x 275mm
Carton Quantity	6pcs/ carton
Carton Weight	9.6kg



HA-552



BMI reading appears immediately after stepping on the scale

- Choose your height before you step onto the scale. Your BMI health range will automatically display together with your weight

Technical Specifications	
Max. Capacity/ Graduation	120kg/ 0.5kg
Product Size (mm)	240 x 280 x 65mm
Product Weight	1.6kg

Packaging Information	
Gift Box Size	244 x 293 x 72mm
Carton Size	305 x 446 x 266mm
Carton Quantity	6pcs/ carton
Carton Weight	11.1kg

Mechanical Scale with Huge Display

Technical Specifications	
Max. Capacity / Graduation	130kg/ 1kg
Product Size (mm)	278 x 292 x 60mm
Product Weight	1.9kg

Packaging Information	
Gift Box Size	283 x 300 x 68mm
Carton Size	312 x 430 x 303mm
Carton Quantity	6pcs/ carton
Carton Weight	13kg

HA-650



HA-680




Classical Mechanical Scale

Technical Specifications	
Max. Capacity / Graduation	130kg/ 1kg
Product Size (mm)	310 x 298 x 64mm
Product Weight	2.1kg
Color Option	White, Black, Yellow

Packaging Information	
Gift Box Size	303 x 318 x 68mm
Carton Size	428 x 350 x 321mm
Carton Quantity	6pcs/ carton
Carton Weight	14.6kg

HA-880

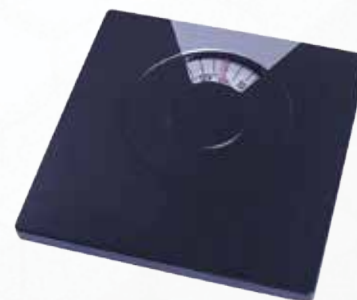
















Mechanical Scale

Technical Specifications	
Max. Capacity / Graduation	130kg/ 1kg
Product Size (mm)	270 x 272 x 50mm
Product Weight	1.7kg
Color Option	White, Black, Yellow

Packaging Information	
Gift Box Size	275 x 280 x 56mm
Carton Size	353 x 292 x 293mm
Carton Quantity	6pcs/ carton
Carton Weight	11.5kg

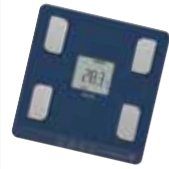Product Comparison for Major Body Composition Monitors

															
	RD-545	RD-953	BC-402	BC-333L	BC-545N	BC-313	BC-G02	BC-541N	BC-730	BC-587	BC-731	BC-722	BC-758	BC-759	BC-770
Max Capacity	200kg	200kg	150kg	150kg	150kg	150kg	150kg	150kg	150kg	200kg	150kg	150kg	150kg	150kg	150kg
Graduation (0-100kg)	0.05kg	0.05kg	0.1kg	0.05kg	0.1kg	0.05kg	0.1kg	0.1kg	0.1kg	0.1kg	0.1kg	0.1kg	0.1kg	0.1kg	0.1kg
Weight	√	√	√	√	√	√	√	√	√	√	√	√	√	√	√
Body Fat %	√*	√	√	√	√*	√	√	√	√	√	√	√	√	√	√
Total Body Water	√	√	√		√	√	√	√	√	√	√	√			
Muscle Mass	√*	√	√	√	√*	√	√	√	√	√	√	√	√	√	√
Muscle Quality	√*	√													
BMR	√	√	√	√	√	√	√	√	√	√	√	√	√	√	
Metabolic Age	√	√	√	√	√	√	√	√	√	√	√	√	√	√	√
Visceral Fat Rating	√	√	√	√	√	√	√	√	√	√	√	√	√	√	√
Physique Rating	√	√	√		√	√		√	√	√	√				
Bone Mass	√	√	√	√	√	√	√	√	√	√	√				
BMI	√	√	√	√	√	√	√					√	√	√	√
Daily Calories Intake															
FitPLUS							√								
Segmental Measurements	√				√										
Auto Recognition	√	√	√	√	√	√	√	√	√			√	√	√	√
No. of User Memories	4 Users	4 Users	5 Users	5 Users	5 Users	5 Users	2 Users	5 Users	5 Users	4 Users	5 Users	5 Users	5 Users	5 Users	5 Users
Product Size/ Weight	348 x 320 x 57mm 2.8kg	298 x 328 x 32mm 2.1kg	217 x 316 x 27mm 1.5kg	274 x 310 x 36mm 1.5kg	320 x 384 x 57mm 2.5kg	217 x 316 x 27mm 1.5kg	226 x 266 x 34mm 0.9kg	285 x 278 x 36mm 1.2kg	216 x 260 x 33mm 880g	300 x 340 x 43mm 2.7kg	210 x 297 x 25mm 1.4kg	215 x 259 x 33mm 1kg	293 x 286 x 25mm 1.7kg	216 x 260 x 33mm 880g	268 x 268 x 35mm 1.1kg
Carton Size/ Weight	381 x 318 x 368mm 14.3kg	265 x 391 x 377mm 14kg	233 x 368 x 278mm 9.5kg	330 x 240 x 338mm 9.1kg	465 x 298 x 398mm 14.4kg	233 x 368 x 278mm 9.5kg	270 x 330 x 320mm 8.3kg	240 x 370 x 365mm 7kg	316 x 288 x 274mm 7.2kg	306 x 410 x 390mm 14kg	223 x 349 x 271mm 8.3kg	245 x 350 x 311mm 7.1kg	228 x 347 x 362mm 11.1kg	316 x 288 x 274mm 7.2kg	302 x 333 x 320mm 8.4kg
Carton Quantity	4 pcs	5 pcs	5 pcs	5 pcs	4 pcs	5 pcs	4 pcs	4 pcs	6 pcs	4 pcs	5 pcs	5 pcs	5 pcs	6 pcs	6 pcs
Tanita Factory Location	Japan	China	China	Japan	Japan	China	China	Japan	China	China	China	China	China	Japan	China

* Segmental Readings Available

TANITA Fitness Gears Your exercise partner!



SR-901



Six Pack / Subcutaneous Fat Analyser

No more Caliper – Analyze your subcutaneous fat level in sec!

- Analyze the thickness of subcutaneous fat around the abdomen, under arm, thigh area etc
- Easy-operation: Moist the skin and place the Analyser, the result will be available in less than 10 secs
- Results will be displayed in both Number and easy-to-read 5-stage Indicator
- The lower the Number, the thinner the subcutaneous fat




Technical Specifications

Product Size (mm)	44 x 26 x 138mm
Product Weight	90g
Batteries	CR2032 x 2

Packaging Information

Gift Box Size	123 x 33 x 161mm
Carton Size	448 x 295 x 225mm
Carton Quantity	24pcs/ carton
Carton Weight	4.8kg

5-stage Indicator:

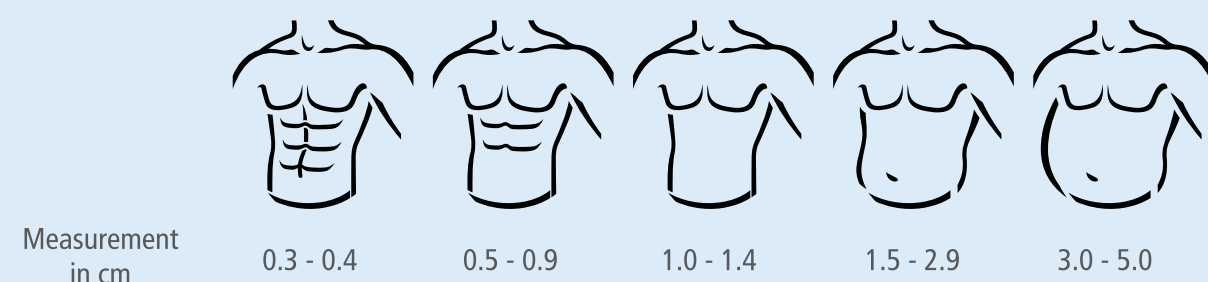

TS-951



Yoga mat

- Non-slippery material with special anti-stretching treatment for best safety
- Light-weight material with convenient carrying bag for handy use
- Washable, hygienic mat material
- Suggested exercise manual attached for quick reference



Not for sale

TS-952



Gym ball

- 65cm large, perfect for 160 – 175cm height
- Can adjust height by altering the amount of air injected
- Special anti-burst treatment for best safety
- Carrying bag, air pump and suggested exercise manual attached



Not for sale

TS-959



Not for sale

Balance Cushion

- Put it on a chair and sit down to train
- Easy-to-use in office
- Suggested exercise manual attached for quick reference



NEW

TANITA Coffee



- A mix of best quality Arabica coffee beans
- Special roasting method to keep the maximum amount of Chlorogenic Acid
- Healthy and tasty



TS-960



Not for sale

Jump Rope with Calorie Counter

- Includes wireless jump rope for using at indoor
- Available for jumps and calories consumed counting



NEW

Elastic Band



Not for sale

- Made of natural Latex, ensuring comfort and durability
- With a supplied storage bag, elastic band can be carried anywhere
- Suggested exercise manual attached for quick reference





The global leader in
precision weighing and BIA body composition equipment

TANITA Corporation

Address: 14-2, 1-Chome, Maeno-cho, Itabashi-ku, Tokyo, JAPAN 174-8630
Tel: +81-3-3968-7048
Fax: +81-3-3968-2661
Email: info@tanita.co.jp
Website: www.tanita.co.jp

TANITA Health Equipment H.K. Ltd

Address: Unit 301-303, 3/F, Wing On Plaza, 62 Mody Road, Tsimshatsui East, Kowloon, HONG KONG
Tel: +852-2834-3917
Fax: +852-2838-8667
Email: asia-pacific@tanita.co.jp
Website: www.tanita.asia

TANITA Europe B.V.

Address: Hoogoorddreef 56E, 1101 BE Amsterdam, THE NETHERLANDS
Tel: +31-20-560-29-70
Fax: +31-20-560-29-88
Email: info@tanita.eu
Website: www.tanita.eu

TANITA Corporation of America, Inc.

Address: 2625 South Clearbrook Drive, Arlington Heights, IL 60005, USA
Tel: +1-847-640-9241
Fax: +1-847-640-9261
Email: 4health@tanita.com
Website: www.tanita.com

TANITA (Shanghai) Trading Co. Ltd.

Address: Room 8005, 887 Huai Hai Zhong Lu, Shanghai, CHINA
Tel: +86-21-6474-6803
Fax: +86-21-6474-7901
Email: info@tanita.com.cn
Website: www.tanita.com.cn

Disclaimer:

TANITA reserved all rights including data integrity of information in this brochure. Shall any parties are interested in the products, please contact TANITA Asia Pacific Regional Sales Office formore information.

TANITA Asia Pacific © Aug 2021 version

Detailed Product Comparison:

Please browse our website

